

## How to do a Twin Bag™ Exchange

### Important things to remember:

- Masks
- Handwashing
- Checking your supplies
- Connection points
- Checking the drained fluid
- Completing your exchange in a clean area with no clutter or drafts


Cũng có sẵn bằng tiếng Việt  
Also available in Vietnamese: VI85-2390

Также есть на русском языке  
Also available in Russian: RU85-2379

同时提供简体中文版  
Also available in Simplified Chinese: CH85-2214




**1**

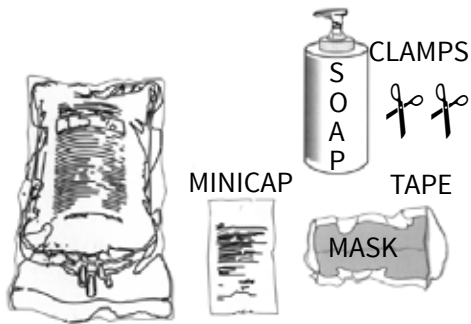


- Mask and wash hands.
- Gather supplies on a clean surface.
- Check Twin Bag™ for: solution strength, bag size, expiry date, and calcium.
- Check MiniCap expiry date.

**2**

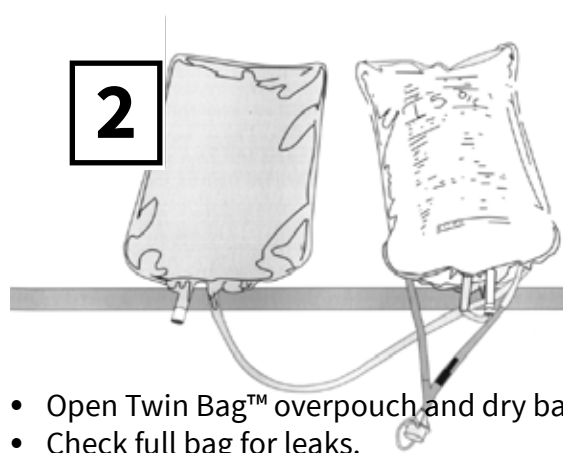


- Open Twin Bag™ overpouch and dry bag.
- Check full bag for leaks.
- Check medication ports.
- Check empty bag.
- Check pull ring, green seal, and Y area.

**3**

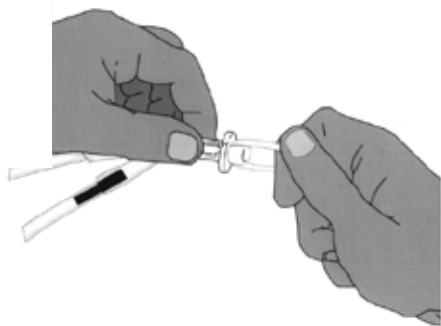


- Tape tubing to the table at the Y.
- Expose transfer set and check that clamp is closed.
- Mask and wash hands with gel.
- Remove the pull ring.

**4**

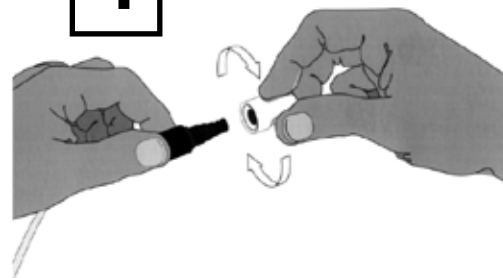


- Remove MiniCap from transfer set.
- Attach Twin Bag™ connector to transfer set right away.

**5**

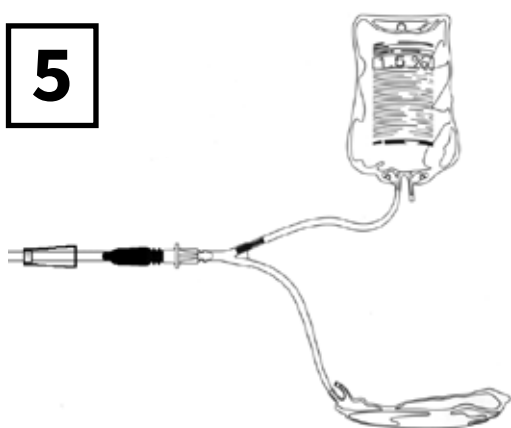


- Hang fill bag and place drain bag down on the pole.

**6**

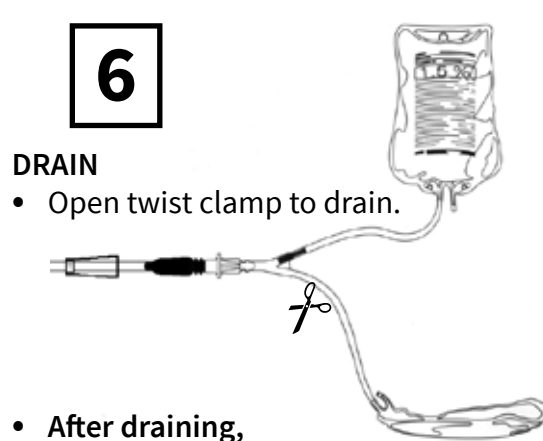


- Open twist clamp to drain.
- After draining, close twist clamp on transfer set.
- Place clamp on drain line.
- Place drain bag on the table to check it.

**7**

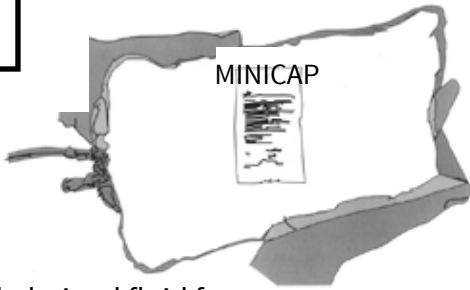

- Check drained fluid for:
  - Clearness: If cloudy, call your Home Dialysis Unit.
  - Fibrin: If fibrin, add heparin as instructed.
  - Blood: If you also have pain, call your Home Dialysis Unit.

**8**

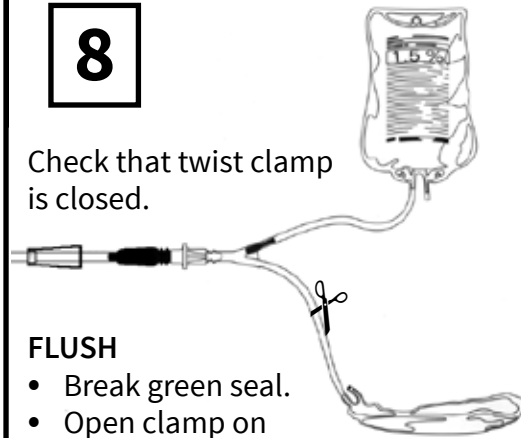

Check that twist clamp is closed.


**FLUSH**

- Break green seal.
- Open clamp on drain line for 3 seconds.
- Close clamp on drain line.

**9**

**FILL**


- Open twist clamp to fill.

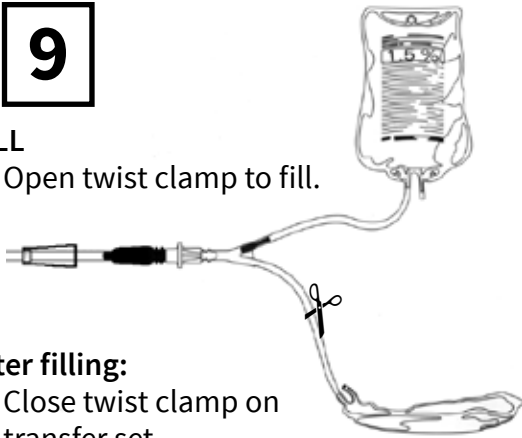

**After filling:**


- Close twist clamp on transfer set.
- Place clamp on fill line.
- Check that all 3 clamps are closed.

**10**

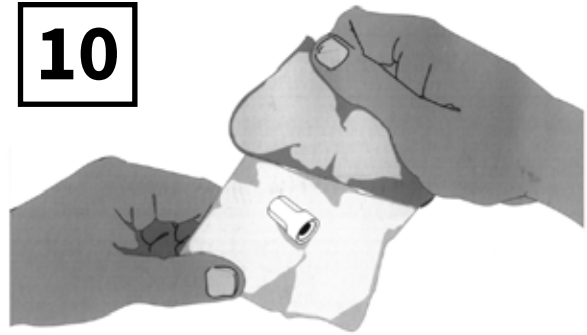


- Mask.
- Wash hands with gel.
- Open MiniCap package.

**11**

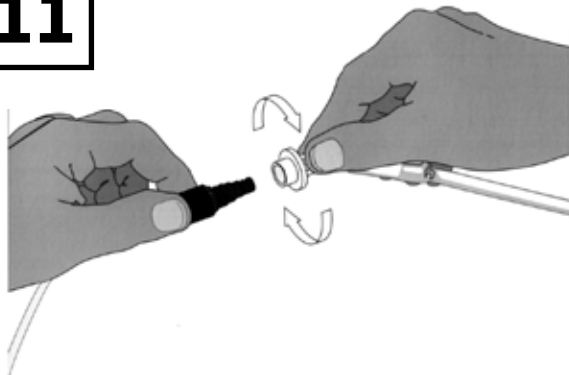

- Unscrew Twin Bag™ connector from transfer set. Keep the transfer set in your hand.

**12**

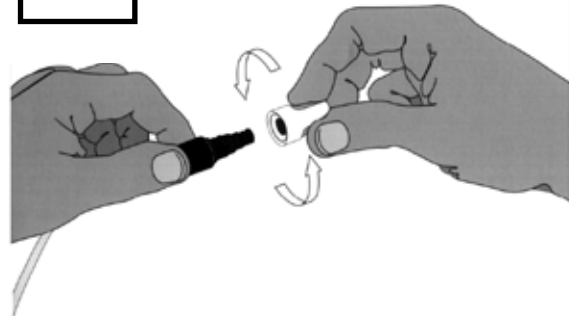

- Apply new MiniCap to transfer set right away.
- Weigh drain bag and record.

# Drain problems - Twin Bag™

<b>Problem/cause:</b>	<b>What to do:</b>
1. Twist clamp	<ul style="list-style-type: none"><li>• Check: is it open?</li><li>• Check under clamp for pinched tubing.</li></ul>
2. Kink in tubing	<ul style="list-style-type: none"><li>• Check around dressing.</li><li>• Check along drain line.</li><li>• Check tape on pole.</li></ul>
3. Fibrin	<ul style="list-style-type: none"><li>• Wrap tubing tightly around pencil and release. (Repeat.)</li></ul>
4. Constipation (not able to poop)	<ul style="list-style-type: none"><li>• Take laxatives as directed by your health care provider.</li></ul>
5. Body position	<ul style="list-style-type: none"><li>• Bend back and forth.</li><li>• Stand up or press on abdomen (stomach area).</li><li>• Cough.</li><li>• Make sure the drain bag is lower than your abdomen.</li></ul>
6. Air lock	<ul style="list-style-type: none"><li>• Wrap tubing tightly around a pencil and release. (Repeat.)</li></ul>
7. Defective bag	<ul style="list-style-type: none"><li>• Change to a new set.</li></ul>

## Fill problems - Twin Bag™

Problem/cause:	What to do:
1. Twist clamp	<ul style="list-style-type: none"> <li>• Check: is it open?</li> <li>• Check under clamp for pinched tubing.</li> </ul>
2. Kink in tubing	<ul style="list-style-type: none"> <li>• Check around dressing.</li> <li>• Check along tubing.</li> </ul>
3. Green seal not broken	<ul style="list-style-type: none"> <li>• Make sure green seal is broken. You should be able to see a space where you cracked the green seal.</li> </ul>
4. Body position	<ul style="list-style-type: none"> <li>• Bend back and forth.</li> <li>• Stand up or press on abdomen.</li> <li>• Cough.</li> </ul>
5. Fibrin	<ul style="list-style-type: none"> <li>• Wrap tubing tightly around a pencil and release. (Repeat.)</li> <li>• Inject heparin in fill bag.</li> <li>• Squeeze bag.</li> <li>• Raise pole.</li> </ul>
6. Constipation	<ul style="list-style-type: none"> <li>• Take laxatives as directed by your health care provider.</li> </ul>
7. Defective bag	<ul style="list-style-type: none"> <li>• Change to a new set.</li> </ul>

## Mechanical problems - Twin Bag™

Problem/cause:	What to do:
1. Transfer set contaminated with twist clamp closed	<ul style="list-style-type: none"> <li>Put MiniCap on for 5 minutes, then remove and put on another new MiniCap.</li> </ul>
2. Did not close twist clamp on transfer set when disconnecting from Twin Bag™	<ul style="list-style-type: none"> <li>Close twist clamp on transfer set, put MiniCap on for 5 minutes, then remove and put on another MiniCap.</li> <li>Call the Home Dialysis Unit <b>RIGHT AWAY. This is a dialysis emergency.</b> There is a nurse on call 24 hours a day, 7 days a week.</li> </ul>
3. Hole in drain bag	<ul style="list-style-type: none"> <li>Close twist clamp on transfer set, then clamp drain line. Put MiniCap on transfer set.</li> <li>Take out new Twin Bag™ system and continue with dialysis exchange.</li> </ul>
4. Hole in fill bag	<ul style="list-style-type: none"> <li>Clamp twist clamp on transfer set.</li> <li>Clamp fill line.</li> <li>Put MiniCap on transfer set.</li> <li>Call the Home Dialysis Unit <b>RIGHT AWAY. This is a dialysis emergency.</b> There is a nurse on call 24 hours a day, 7 days a week.</li> </ul>
5. Your tubing separates from your catheter: <ul style="list-style-type: none"> <li>› Connection points are not tight</li> <li>› There is pulling on catheter and tubing</li> </ul>	<ul style="list-style-type: none"> <li>Use emergency kit. Put catheter clamp on your catheter next to the exit site.</li> <li>Put cap on end of your catheter.</li> <li>Call the Home Dialysis Unit <b>RIGHT AWAY. This is a dialysis emergency.</b> There is a nurse on call 24 hours a day, 7 days a week.</li> </ul>
6. Needle goes through port of new bag	<ul style="list-style-type: none"> <li>Change to new Twin Bag™ set and new syringe and needle.</li> <li>Draw up new medication.</li> </ul>
7. New fluid goes in drained bag	<ul style="list-style-type: none"> <li>Do another bag exchange.</li> </ul>
8. Contaminated Twin Bag™ connector	<ul style="list-style-type: none"> <li>Throw away the set.</li> <li>Check and prepare a new set.</li> </ul>



This pamphlet is for educational purposes only. It is not intended to replace the advice or professional judgment of a health care provider. The information may not apply to all situations. If you have any questions, please ask your health care provider.

Find this pamphlet and all our patient resources here:  
<https://library.nshealth.ca/Patients-Guides>

Connect with a registered nurse in Nova Scotia any time:  
Call 811 or visit: <https://811.novascotia.ca>

*Prepared by:* Renal Program  
*Illustrations used with permission from:* Baxter Health Care ©  
*Designed by:* Nova Scotia Health Library Services

WQ85-0633 © January 2025 Nova Scotia Health Authority  
To be reviewed January 2028 or sooner, if needed.

